

CRISPOR Manual

Maximilian Haeussler, Jean-Paul Concordet

Download this page as a PDF or Epub. Comments or questions? Email us at crispor@tefor.net.

What is CRISPOR?


CRISPOR is a website that helps select and express CRISPR guide sequences. In its default mode, the user pastes an input DNA sequence and chooses the genome. CRISPOR then lists the guides in the input sequence, adding information about them found in databases and algorithms, including genome variants, predicted off-target and on-target scores. For each guide sequence, various primers are designed, e.g. to amplify the target, produce guide RNAs by in vitro transcription after annealing of overlapping primers or for cloning into AddGene plasmids.


A second, optional mode, called CRISPOR Batch, is available for users who want to use pre-selected guide RNAs for gene inactivation experiments in mouse or human cells. It accepts one or multiple gene identifiers and returns a number of pre-selected guides from various genome-wide libraries as well as primers with overhangs for cloning into guide RNA expression plasmids and with sequence information preformatted for easy ordering from oligo-array manufacturers.

Additional assistants in CRISPOR design primers for all predicted off-targets of a guide and create oligonucleotide pools for saturating mutagenesis of the input sequence.

Below, we provide a step by step manual of CRISPOR and instructions on how to use the assistants. Important information is also directly available on the web pages by hovering the mouse over the small info bubbles shown in orange. Usually, orange colour indicates links which can be clicked for further information or more detailed CRISPOR results.

Input


CRISPOR is a program that helps design, evaluate and clone guide sequences for the CRISPR/Cas9 system. 

Step 1
Planning a lentiviral gene knockout screen? Use [CRISPOR Batch](#)
Sequence name (optional):
Enter a single genomic sequence, < 2000 bp, typically an exon 
[Clear Box](#) - [Reset to default](#)

```
cttctcttgtccccaattggggcgccggcccccctggggcctaaaggactggcgccggaagtgccacgg  
ggggggggaactcggcgcacagcgccggcctattcttcggcgtacacgatgatggaatcgcgccgc  
tcctctctcgacacggcctccgscatgtgcaaggccggcttcgcgggccgacgatgcccccc  
cgggccgtcttccccctcatcgtggggcgcc
```

Text case is preserved, e.g. you can mark ATGs with lowercase.
Instead of a sequence, you can paste a chromosome range, e.g. chr1:11908-12378

Step 2
Select a genome
Homo sapiens - Human - UCSC Feb. 2009 (GRCh37/hg19) + SNPs: 1000Genomes, ExaC
Note: pre-calculated exonic guides for this species are on the [UCSC Genome Browser](#).
We have 246 genomes, but not the one you need? Send its FASTA/GFF URL to [CRISPOR support](#)

Step 3 
Select a Protospacer Adjacent Motif (PAM)
20bp-NGG - Sp Cas9, SpCas9-HF1, eSpCas9 1.1

SUBMIT

Version 4.3, Feedback: [email](#) - [Downloads/local installation](#)

Figure 1: The sequence input page of crispor.org

On the first page of CRISPOR.org, you have to enter only three pieces of information:

1. **Input Sequence:** paste your target sequence into the input sequence box. No special sequence format is required (FASTA accepted). Characters other than A, C, G, T and N will be automatically removed from the pasted sequence. A warning will be displayed during processing if characters were removed.
 - the input sequence is usually genomic DNA sequence but any input sequence can be submitted
 - your input sequence should usually be contained in the selected genome, but it may not be present, e.g. if you are designing guides against a transgene, like GFP. A warning will appear in this case. This warning will also appear if you have accidentally used cDNA sequence as input. Using cDNA is not recommended since guide sequences that overlap exon-exon boundaries will not have bona fide targets in the genome.
 - case of the input sequence will be retained in the output. As such, you can mark any sequence of interest (exons, ATG or stop codons, coding sequence of catalytic sites, . . .) using upper/lower case in order to easily visualize where the different guide sequences are positioned relative to these features in the output.
 - if the input sequence contains N-characters (“aNy” nucleotide), no guides will go over these characters, so you can mark positions that you want to exclude from the design with Ns, e.g. to avoid single nucleotide polymorphisms (SNPs also known as SNVs, single nucleotide variants).
 - you can give your input sequence an optional name. It will be

shown on the output page and appended to all oligonucleotides that you download later.

2. Genome: Select your genome of interest from the list. You can scroll through and select from the drop-down or type and search the names of the species, using latin or English common names.
 - a default genome will be automatically set to the last genome that you selected.
 - we have imported numerous genomes already. Some species are available multiple times, e.g mouse and human, because assemblies from different years are available and some loci are only part of certain assembly releases. Also, the annotation with variants (SNPs and short indels, see below) is only available for certain assemblies. This is shown as part of the genome list, e.g. the 1000 Genomes variant annotation is only available for the human genome assembly called “hg19” (aka GrCh37).
 - if your genome is not on the list, please contact us and send us a link to the fasta file and ideally also GFF gene annotations and the common and scientific names, e.g. “zebrafish” and “Danio rerio”. If your genome is in NCBI RefSeq or on the UCSC or Ensembl browsers, please send us the NCBI assembly accession ID or a link to the UCSC or Ensembl page. We have added more than 150 genomes since the publication of the CRISPOR paper.
3. Protospacer Adjacent Motif (PAM): For most current applications of the CRISPR-Cas system, *Streptococcus pyogenes* Cas9 nuclease is used and the corresponding PAM is NGG, the default. However, you can choose other enzymes and corresponding PAMs from the dropdown box.

Output 1: Annotated input sequence

The main output of CRISPOR is a page that shows the annotated input sequence at the top and the list of possible guides in the input sequence at the bottom.

The input sequence is shown first. Underneath the sequence, all PAM (Protospacer adjacent motif) sites are highlighted. Most labs use Sp-Cas9, its PAM is NGG. Sp-Cas9 usually cuts three or four basepairs 5' of the PAM site. The three bases 5' of the PAM are marked with dashes (“-”) to help users find the cleavage position. PAM sites can also be on the reverse strand of the input sequence. For Sp-Cas9, these correspond to CCN motifs on the input sequence and the dashes are shown on the right side of the PAM in

[Homo sapiens \(hg 19\)](#), [chr7:5569177-5569415](#), [reverse genomic strand](#) [Link to Browsers](#)

Found 52 possible guide sequences in input (239 bp). Click on a PAM NGG match to show its 20 bp guide sequence.

Shown below are the PAM site and the expected cleavage position located -3bp 5' of the PAM site.

Colors **green**, **yellow** and **red** indicate high, medium and low specificity of the PAM's guide sequence in the genome.

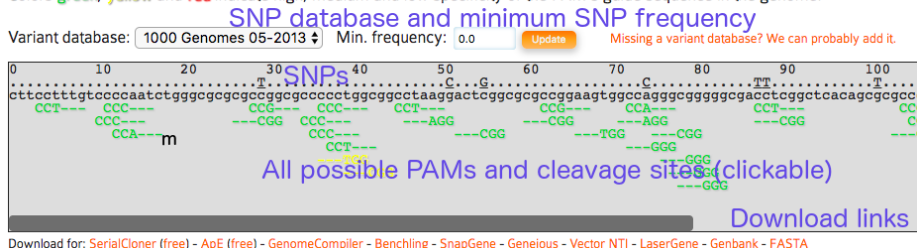


Figure 2: The annotated sequence

this case. PAMs are clickable and link to the corresponding target row in the table of guides below (see next section).

- you can send the URL of the output page to collaborators or bookmark it. We have not yet deleted any results over the last two years and try to keep the output for at least a year. For long-term archival, consider downloading the page and/or the Excel spreadsheets.
- the position of the input sequence in the selected genome is shown at the top of the main CRISPOR results page. A click onto this position will open a new window with a genome browser, usually UCSC or Ensembl, depending on the source of the genome. If the input sequence was not found in the selected genome, a warning is displayed (see note in section here above).
- for certain genomes, the input sequence can be annotated with variants, mostly SNPs. The variants are shown above the input sequence. You can hover with your mouse over them, to show details about the variant, usually the nucleotide change (T>G) and the frequency. You can change the variant database and you can also set a minimum frequency for variants shown on the page. Any variant that has a frequency below the threshold will not be displayed.
- very few genomes (human and mouse) have variants in our database at the moment. If you need a particular one, contact us by email and send us the URL of the database or the VCF file.

Download as Excel tables: [Guides](#) / [Off-targets](#) / [Saturating mutagenesis assistant](#)

Position/ Strand	Guide Sequence + PAM + Restriction Enzymes + Variants <input type="checkbox"/> Only G- <input type="checkbox"/> Only GG- <input type="checkbox"/> Only A-	Specificity Score	Predicted Efficiency Show all scores Doench -16 Mor-Mateos		Out-of-Frame score Click on score to show micro- homology	Off-targets for 0-1-2-3-4 mismatches + next to PAM	Genome Browser links to matches sorted by CFD off-target score <input type="checkbox"/> exons only <input type="checkbox"/> chr7 only
46 / rev	TTCGGCGCGCGAGTCCTT AGG Enzymes: BshFI, BseZ1I, TauI, BstDEI, SstI, Fsp4II Cloning / PCR primers	97	49	57	67	0-0-0-0-15 0-0-0-0-1 15 off-targets	4:intergenic:SYNJ2BP/SYNJ2BP-COX16- RP11-486O13.4 4:intergenic:AC013275.2-SCTR 4:intron:CR1L show all...
12 / rev	CGCGCGCGCGCCAGATT GGG High GC content Enzymes: BstFI, BstI Cloning / PCR primers	96	21	58	71	0-0-0-3-28 0-0-0-0-0 31 off-targets	4:exon:LMP1 4:intergenic:MIR4291-RP11-53B5.1 4:intron:DAZAP1 show all...
30 / fw	TCCCAATCTGGCGCGCGC CGG Enzymes: NlaIV, AclI, BfoI, MspI, BssHII,	93	37	36	65	0-0-0-0-55 0-0-0-0-1	4:intron:EPS15L1 4:exon:USP32 4:exon:SBF1

Figure 3: The guide list

Output 2: Guide list

Shown below the input sequence are the guide target sequences, one per PAM. For spCas9, the PAM is NGG and the targets are 20bp long. Each 20bp target sequence in the input sequence is aligned against the whole genome allowing at most four mismatches and the results are summarized as a table. The table has the columns described below. To sort by one of them, click its name in the first row.

By default, the table rows are sorted by specificity score. To sort by values of any other column, click the name of the column in the first row.

Column 1 - guide name: this is the position of the PAM on the input sequence and the strand, e.g. "13+". You can sort by position by clicking the header "Position/Strand" in the table.

Column 2 - guide sequence: the sequence of the guide target and the PAM and also one of the most important features of CRISPOR, the link to its "PCR and cloning primers" (see the Primers section below).

In addition, depending on the genome and guide, additional data is displayed in the guide sequence column:

- variants are indicated underneath the target sequence and PAM, if variants are available for this genome.
- High-GC-content (>80%) guides are flagged, as are low-GC content (<20%) targets. Various studies have reported that both cases lead to low target cleavage efficiency.
- Restriction enzyme sites that overlap the three base pairs 5' of the PAM site. These will likely be disrupted during DNA repair of the double-strand break made by the Cas9 nuclease and may therefore be

convenient to use for screening of mutant sequences amplified by PCR.

- Also shown in this column are restriction enzymes that overlap the three basepairs 5' of the PAM site. Finally, one of the most important features of CRISPOR is available here, under the link "PCR primers" .

This column allows to filter the guides by the first nucleotide, as some RNA expression promoters can only transcribe guides that start with certain nucleotides. Some labs only use targets that start with G- for the U6 promoter, A- for the U3 promoter or G- for the T7 promoter. In this way, no further sequence changes are necessary. In our lab, we do not constrain guides in this way, but simply prefix them with the required nucleotide as a 21st basepair (see the Primers section below).

Column 3 - specificity score: this score is a prediction of how much an RNA guide sequence for this target may lead to off-target cleavage somewhere else in the genome. The score ranges from 0-100 with 100 being the best, meaning the search could not find a single sequence in the genome that differs from the target at up to four positions. This score uses the formula from the MIT Crispr Website (Hsu guide score) but with a better and more sensitive search engine. We think that good guides should have a specificity score of at least 50, based on the data from whole-genome off-target assays, see Figure 3a in the CRISPOR paper. The color of the guide (red/yellow/green) is based on the specificity score (>50 = green, >30 = yellow). You should avoid guides with very low scores, unless you can validate the off-targets with special assays or you can cross or inbreed the animals until you are sure that no off-targets are left. SaCas9 now has its own specificity score provided by Josh Tycko, so this column will appear for SaCas9, but there is no such score for Cpf1 yet.

Column 4 - efficiency scores: the efficiency score is a prediction of how well this target may be cut by its RNA guide sequence. It ranges from 0-100 with 100 being the best. We show two scores here: the scoring method from the Doench 2016 paper (aka "sgRNA Designer") or the one by Moreno-Mateos 2016 ("CrisprScan"). In our study, we found that the Doench 2016 score is the best score for guides expressed in the cells from a U6 promoter, see Figures 4 and 5 in our CRISPOR paper. The Moreno-Mateos2016 score was better in our study when the guide was expressed in-vitro with a T7 promoter. While there is link in the header of this column to show various other scores in the table, we do not recommend the other algorithms anymore, based on the results from our paper. Note that the Doench 2016 scores slightly changed in April 2018, when we updated to the version called 'Azimuth', which is also used by the Broad Institute's website, but the changes were on the order of 1-2%, so it did not make any difference in practice.

Cpf1 and saCas9 have special efficiency scores, the spCas9 scores are not

applicable for them, and the models “deepCpf1” for Cpf1 and “Najm et al 2018” for saCas9 will be shown here.

Also, the predictive power of these scores is not great in general, with correlations of around 0.4 against the assay results. If you have no choice and need to pick certain guides, you may choose to ignore efficiency scores. If you have a choice between many specific guides, however, or you are screening many guides, you may want to prefer guides with high efficiency scores. For a more detailed discussion of this topic, see our Notes on on-target cleavage below.

Column 5 - out-of-frame score : this score (0-100) is a prediction how likely a guide is to lead to out-of-frame deletions. This is relevant if you are doing gene knockouts with a single guide. Gene knockouts with single guides work because repair after DNA cutting is error-prone and small deletions are introduced 5' of the PAM. It has been observed by Bae et al. that this repair does not lead to a random distribution of small deletions, but that due to microhomology around the cut site, certain deletions are favored, depending on their flanking DNA sequences. You can click on this score to show what the predicted deletions are, in order of their predicted microhomology. We know of at least one case where it was impossible to obtain a gene knock-out with a specific guide because all deletions seen after sequencing were always in-frame (T. Momose, unpublished data). The higher the out-of-frame score, the more deletions have a length that is not a multiple of three, see Bae et al..

Column 6 - off-target mismatch counts: the number of possible off-targets in the genome, for each number of mismatches. This is a summary of the whole-genome search for sequences similar to the guide target sequence. It is best explained by an example: a description “0 - 1 - 2 - 9 - 28” means that the target matches 0 locations in the genome with no mismatch, 1 location in the genome with 1 mismatch, 2 locations with 2 mismatches, 9 with 3 and 28 locations with 4 mismatches. The smaller numbers in grey below use the same scheme, but only for locations with no mismatch in the 12 bp close to the PAM, the “seed” region. Early reports based on in-vitro cleavage assays suggested that off-targets with a mismatches in the seed region are very inefficiently cut. If the grey numbers are “0 - 0 - 1 - 7 - 2” this means that there is one location in the genome where the target matches with three mismatches and all three mismatches are outside the seed region, 7 locations with 4 mismatches outside the seed region and 2 locations with four mismatches outside the seed region. The total number of off-targets is shown in this column, too. For detailed discussion of the topic of off-target searching, see our Notes on off-target cleavage below.

Column 7 - off-targets: here CRISPOR lists the locations of all possible off-

targets with up to four mismatches, annotated with additional information: - genomic position and an annotation whether they fall into an exon, intron or between genes and the closest gene. - you can hover with your mouse over the position to show an alignment of the guide sequence with the off-target and two off-target scores, CFD (Cutting frequency determination) and MIT.

- by default only the three most likely off-targets are shown, click on “show all” to see more.
- when the “show all” is active, a link to “Off-target primers” is visible. This is a tool for bulk primer design, see Off-target primers below.
- it is possible to shown only off-targets in exons or off-targets on the same chromosome as the target. Exonic off-targets are more likely to have a functional impact. When working with experimental organisms, off-targets lying on a different chromosome than the target will not cosegregate and will be diluted out during breeding.
- All off-targets are sorted by an “off-target score” which tries to predict what the most likely off-targets are. The score is based on where the mismatches are. Usually, a mismatch close to PAM means that the off-target is less likely. See our Notes on off-target cleavage below for more details.
- For some genomes, there is no exon filter. This is the case when we could not find a valid gene model (GFF) file for the genome. Do not hesitate to contact us by email if you have a gene model file or want us to add one.
- When searching for off-targets, CRISPOR is somewhat more tolerant and allows not only NGG, but also NGA and NAG, as these have been found to be quite common in the Tsai et al study, up to 10% of off-target events had non-NGG PAMs. Other enzymes have similar flexibility, GAW for NGN and NKG PAMs, NGG for NGA PAMs, and NNGRRN for NNGRRR PAMs. The tool-tip in the off-target column always shows the currently allowed non-canonical PAMs tolerated in the off-target search phase.

Output 3: Primers

The third main page of CRISPOR is the list of primers related to a selected guide. You reach it by clicking “Cloning / PCR primers” for a particular guide in the table.

The top part of the page (“Cloning or expression of guide RNA”) shows expression primers. There are many possible ways for guide expression supported by CRISPOR. The choice depends on your organism (mice/zebrafish vs. cell cultures) and preferred cloning method. We tried to link to useful protocols

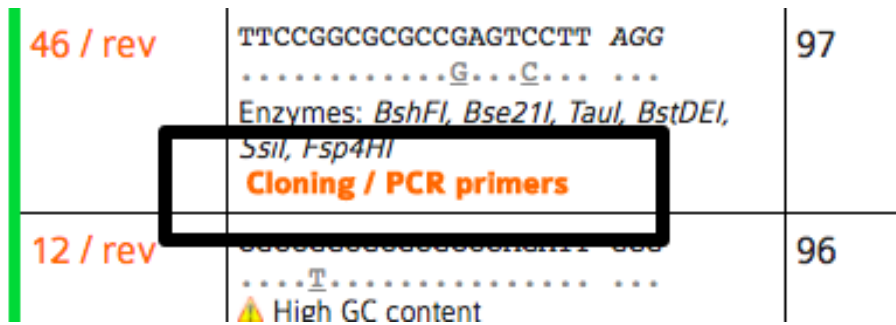


Figure 4: Link to primer designer

Guide sequence: TTCCGGCGCGCCGAGTCCTT AGG

Contents:

- Cloning or expression of guide RNA
 - T7 *in vitro* expression from a plasmid
 - T7 *in vitro* expression from overlapping oligonucleotides
 - U6 expression from an Addgene plasmid
 - Direct PCR for *C. intestinalis*
 - Lentiviral vectors: Cloning with Gibson assembly
 - Summary of main cloning/expression primers
- PCR to amplify the on-target site
- Restriction sites for PCR validation
- PCR to amplify off-target sites
- Saturating mutagenesis using all guides

Figure 5: Contents of the primer designer page

on this page, too. If this page is missing support for your favorite system or protocol, do not hesitate to contact us. For quicker copy/pasting into Excel, this section contains a table with a copy of all primers at the end (“Summary of main cloning/expression primers”).

The second part of this page (“PCR to amplify the on-target site”) shows one pair of primers to amplify a fragment of the genome around the target site, to validate cleavage, e.g. with the T7 endo assay, by analyzing your Sanger traces with TIDE or ICE, NGS sequencing analyzed with Crispresso or by restriction digest (see below). You can change the size of the fragment and the T_m of the primers here to optimize the PCR design to the specific assay.

The third part of this page (“Restriction sites for PCR validation”) lists restriction enzymes that overlap the position 3bp 5′ of the PAM. If a PCR product of the primers above is not cut with the enzyme anymore, it is very likely that the genome had been changed by Cas9 successfully, usually with a short deletion. This is quicker in the lab than sequencing or T7 endonuclease assays and can be very useful for routine screening during breeding of mutant animals.

To make sure that your restriction site is unique enough in the PCR fragment to be visible on a gel, the sites are shown on the PCR fragment sequence shown below. The sizes of the individual sub-fragments are shown, too.

Finally, there are links at the bottom of the page to design primers for all off-targets of this guide (see next section) and also a link to download all guides in the input sequence for saturating mutagenesis, e.g. of non-coding regions (see below).

CRISPOR off-target primer assistant

High-throughput sequencing is the most accurate way of determining off-target cleavage. By ordering oligos and running arrayed PCRs e.g. on 96-well plates, you can validate up to several hundred off-targets relatively quickly. The off-target primer designer of CRISPOR outputs not only one pair of primers for each off-target with the right overhangs for Illumina adapter addition, but also creates an input file for CRISPRESSO. CRISPRESSO (<http://crispresso.rocks>) is a software package for the analysis of FASTQ files from sequencing, it can e.g. plot histograms with the distribution of indels 5′ of the PAM site and calculate p-Values for them.

For off-target amplification, we recommend the protocol by Matt Canver where two PCRs are run: one PCR to amplify the potential off-target, then a

second PCR to extend the handles with Illumina barcodes. You can download the protocol here.

Predicted Efficiency Show all scores		Out-of-Frame score Click on score to show micro-homology	Off-targets for 0-1-2-3-4 mismatches + next to PAM Show all scores	Genome Browser links to matches sorted by CFD off-target score <input type="checkbox"/> exons only <input type="checkbox"/> chr7 only
Doench '16	Mor-Mateos	67	0-0-0-0-15 0-0-0-0-1 15 off-targets	4:intergenic:SYNJ2BP/SYNJ2BP-COX16-RP11-486O13.4 4:intergenic:AC013275.2-SCTR 4:intron:CR1L 4:intergenic:W12-81516E3.1-RPL35P8 4:intergenic:DNAH17-CTD-2357A8.3 4:intergenic:RP11-795F19.5-RP11-795F19.5/RBFADN 4:intergenic:EMX2-CTA-109P11.1 4:exon:C3orf58 4:intergenic:PAIP2B-ZNF638 4:intron:CHID1 4:intron:PLCXD1 4:intergenic:Gap-RNU6-1334P 4:exon:CILP2 4:exon:ATG2B/GSKIP 4:intron:ZADH2 show less...

Figure 6: Link to off-target primer designer

You reach the off-target primer designer by first clicking “show all” in the list of off-targets, then following the link “Off-target primers”. At the top of the page, a table with one pair of primers for every predicted off-target is shown.

Name	Primer Sequence	Tm	CFD Score
ontarget_mm0_exon_ACTB/AC006483.1_chr7_5569347_F	TCGTCGGCAGCGTC GTGTGAGCTGCGAGAATAGC	60.1	1.00
ontarget_mm0_exon_ACTB/AC006483.1_chr7_5569347_R	GTCTCGTGGGCTCGG CGGCTTCCTTTGTCCCAAT	60.6	1.00
mm4_intergenic_SYNJ2BP/SYNJ2BP-COX16 RP11-486O13.4_chr14_70883812_F	TCGTCGGCAGCGTC CACACGGGTTTCGGTTTCAG	60.0	0.49
mm4_intergenic_SYNJ2BP/SYNJ2BP-COX16 RP11-486O13.4_chr14_70883812_R	GTCTCGTGGGCTCGG ACTGCAATCCATTGGCGGTA	60.0	0.49
mm4_intron_CR1L_chr1_207844275_F	TCGTCGGCAGCGTC CACTAGAAATTTTGTGCTGGAGT	58.8	0.21
mm4_intron_CR1L_chr1_207844275_R	GTCTCGTGGGCTCGG TGCATACACCTTGGGAGCAT	59.3	0.21
mm4_intergenic_W12-81516E3.1 RPL35P8_chr22_49591792_F	TCGTCGGCAGCGTC GGCAACCTCCCTTGTCTGG	60.3	0.19
mm4_intergenic_W12-81516E3.1 RPL35P8_chr22_49591792_R	GTCTCGTGGGCTCGG TCACTCGGGAGGAGTGAG	60.0	0.19

Figure 7: List of off-target primers

In the table of primers output by CRISPOR, Illumina Nextera Handle sequences have been added and are highlighted in bold. Primers for the on-target sequence have been added for convenience. The table is sorted by the CFD off-target score. Note that sites with very low CFD scores < 0.02 are unlikely to be cleaved, see Figure 2 of our CRISPOR paper, but they are still shown on this table. All primer names are prefixed by the number of mismatches, e.g. “mm4” is an off-target with four mismatches. Naturally,

each off-target has two primers, indicated by the common suffixes “_F” and “_R”.

Name	Primer Sequence	Tm	CFD Score
ontarget_mm0_exon_ACTB/AC006483.1_chr7_5569347_F	TCGTCGGCAGCGTCTGTTGAGCTGCGAGAAATAGC	60.1	1.00
ontarget_mm0_exon_ACTB/AC006483.1_chr7_5569347_R	GTCTCGTGGGCTCGGCGGCTTCCTTTGTCCCAAT	60.6	1.00
mm4_intergenic_SYNJ2BP/SYNJ2BP-COX16 RP11-486O13.4_chr14_70883812_F	TCGTCGGCAGCGTCCACAGCGGTTTCGGTTTCAG	60.0	0.49
mm4_intergenic_SYNJ2BP/SYNJ2BP-COX16 RP11-486O13.4_chr14_70883812_R	GTCTCGTGGGCTCGGACTGCAATCCATTGGCGGTA	60.0	0.49
mm4_intron_CR1L_chr1_207844275_F	TCGTCGGCAGCGTCACTAGAAATTTTGTGCTGGAGT	58.8	0.21
mm4_intron_CR1L_chr1_207844275_R	GTCTCGTGGGCTCGGTGCATACACCTTGGGAGCAT	59.3	0.21
mm4_intergenic_WI2-81516E3.1 RPL35P8_chr22_49591792_F	TCGTCGGCAGCGTCCGGCAACCTCCCTGTCTGG	60.3	0.19
mm4_intergenic_WI2-81516E3.1 RPL35P8_chr22_49591792_R	GTCTCGTGGGCTCGGTCACTCTCGGAGGAGTGAG	60.0	0.19

Figure 8: List of off-target primers

The primer designer also outputs a table of the amplicons with the putative off-target sequence highlighted in bold.

Off-target amplicon sequences

Primers underlined, off-targets in bold.

ontarget_mm0_exon_ACTB/AC006483.1_chr7_5569347	<u>TGGT</u> GAGCTGCGAGAAATAGC CGGGCGCGCTGTGAGCCGAGTCCGCCCCCTGGCCACT TCGGCGCGCCGAGTCTCT ANNNNNNNNNNNNNNNNNNNNNNNNNNNNNNNAGATTGGGGACAAAGGAAGCCG
mm4_intergenic_SYNJ2BP/SYNJ2BP-COX16 RP11-486O13.4_chr14_70883812	<u>CACAGCGGTTTCGGTTTC</u> CAGCAGCCTCGAGACCCGGAAAGGAAGCCCGAAGACTCGGAGCACTGGGAA TGCTGGGAAT CGTGGTTCT TACCGCAATGGATTGCAGT
mm4_intron_CR1L_chr1_207844275	<u>ACTTAGAAATTTTGTG</u> TGCTGGAGTAATTTCAATTTTAAATGCTTATGGATGCTCCGGTGTGCACAGTCTTTGGAAT TAAAGAAATGTAAGAGAAAGAGCT TATGCTCCCAAGGTGTATGCA

Figure 9: List of off-target amplicons

Finally, at the end of the page, you can download CRISPRESSO input files and the exact command line to analyze your FASTQ sequencing result file.

Input file for Crispresso

Crispresso, written by Luca Pinello, is a software package to quantify the Cas9-induced mutations on off- or on-targets.

[Click here](#) to download an amplicon input file for Crispresso. For each off-target, it includes the off-target name, its PCR amplicon and the guide sequence. Keep a copy of this file.

After sequencing, run CRISPRessoPooled. The tool will map the reads to the amplicons and analyse the mutations:
 CRISPRessoPooled -r1 Reads1.fastq.gz -r2 Reads2.fastq.gz -f crisporAmplicons_3TeBYMvGLsICUasefw16.txt --name MY_EXPERIMENT

Figure 10: Link to download Crispresso files and sample command line

CRISPOR saturating mutagenesis assistant

In a saturating mutagenesis experiment, a target region of the genome is modified with many guides, to create as many DNA edits as possible. Cells can be transfected with a pool of oligonucleotides (“pooled”) or one guide

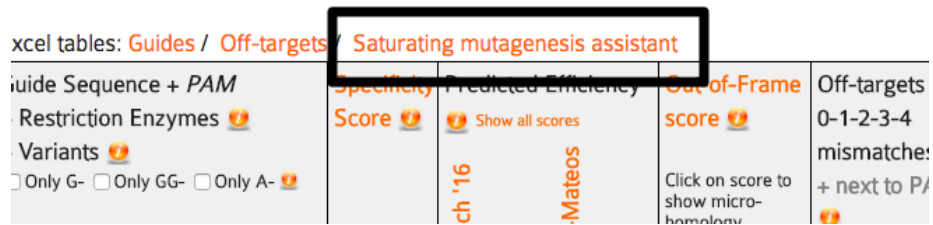


Figure 11: Link to Saturation Mutagenesis Assistant from the list of guides

per cell sample on plates (“arrayed”) and are then phenotyped. This assay is usually used for cis-regulatory regions, to find the relevant sequences for a phenotype. Pools with up to several thousand guides can be ordered from custom oligonucleotide array manufacturers, containing subpools of guides that can induce thousands of deletions into a target sequence.

The CRISPOR Saturating mutagenesis assistant makes this much easier by generating tables with all relevant guides in a target sequence, prefixed with the right adapter sequences, input files for CRISPRESSO to quantify the guides (see off-target primers) and sequencing primers to validate the modifications.

For step-by-step details on how to conduct this assay and use this assistant, please refer to our pre-print, Canver et al. 2017, to appear in Nature Protocols in 2018.

CRISPOR Batch

Pooled CRISPR screens can be used to knock-out thousands of genes with a single experiment. CRISPOR Batch makes the design of an oligonucleotide pool very easy: paste a list of genes, then download a table of oligonucleotides. We currently only offer this for human and mouse, let us know if you need another organism. We also do not design the oligonucleotide libraries but instead offer existing and tested sequences from published libraries. To reach CRISPOR batch, follow the link to it from the main page.

Select the CRISPR library, the number of negative controls, the number of guides you need per gene, then paste a list of gene symbols, Entrez Gene IDs or Refseq IDs, one per line into the box below.

The output page contains a table with one row per guide, for up to the number of guides per gene that you selected.

A note on money: pooled screens may seem like an expensive assay, but at

Step 1

Planning a lentiviral gene knockout screen? Use [CRISPOR Batch](#)

Sequence name (optional):

Enter a single genomic sequence, < 2000 bp, typically an exon

[Clear Box](#) - [Reset to default](#)

```
cttcctttgtccccaatctggggcgcgccggcgc
```

Figure 12: Link to CRISPOR Batch from the main page

Lentiviral screen library:

Number of guides per gene:

Number of non-targeting control guides (max: 1000):

Subpool barcode:

Enter a list of gene symbols, Entrez Gene IDs or Refseq IDs, one per line (case-insensitive):
 You may need to [convert old symbols](#). Type 'all' below to get all guides in the library and no gene filtering

```
PITX2
MTOR
TP53
ABO
3861
NM_134261
```

Figure 13: CRISPOR Batch Input

Library: Human, Brunello, Doench Nat Bio 2016
Number of guides per gene: 3
Number of non-targeting controls: 10 (all controls are from the GeCKOV2 library)
Subpool Barcode: Subpool 1: CGGGTTCCGT/GCTTAGAATAGAA
Subpool Prefix: CGGGTTCCGT
pLentiGuidePre primer: GGAAAGG
pLentiGuidePre: ACGAAACACCG
pLentiGuidePost: GTTTTAGAGCTAGAAATA
pLentiGuidePost primer: GCAAGTTAAAATAAGGC
Subpool Suffix: GCTTAGAATAGAA
Oligonucl. structure: Subpool Prefix + pLentiGuidePre Primer + pLentiGuidePre + sgRNA + pLentiGuide Post + pLentiGuidePost Primer + Subpool Suffix

For details on these sequences, see our [protocol](#)

Input gene identifiers that were not found: NM_134261

[Download table](#)

ID of guide	Target Entrez ID	Target Refseq ID	Guide RNA (click to show in CRISPOR)	Full oligonucleotide including guide RNA
PITX2_1	5308	NM_000325.5	CCCTACGACGACATGTACCC	CGGGTTCCGTGGAAAGGACGAAACACCGCCCTACGACGACATGTACCCGTTTTAGAGCTAGAAATAGC
PITX2_2	5308	NM_000325.5	CACAGCGATTTCTTCGCGTG	CGGGTTCCGTGGAAAGGACGAAACACCGCACAGCGATTTCTTCGCGTGGTTTTAGAGCTAGAAATAGC
PITX2_3	5308	NM_000325.5	AGCAGCGGGGACCAAGCTGT	CGGGTTCCGTGGAAAGGACGAAACACCGAGCAGCGGGGACCAAGCTGTGTTTTAGAGCTAGAAATAGC
MTOR_1	2475	NM_004958.3	TCAGGAAATGATCCGCACAG	CGGGTTCCGTGGAAAGGACGAAACACCGTCAGGAAATGATCCGCACAGTTTTAGAGCTAGAAATAGC

Figure 14: CRISPOR Batch Output

current prices you can get twelve thousand oligos for < 2000\$ and the pool can easily be split into smaller subpools by PCR. In addition, there is a lot of competition in this field. Our collaborator Matt Canver uses CustomArray which says they will beat any quote. If you order through Twist Biosciences, mentioning CRISPOR will get you 35% off your order (we do not receive kick-backs from them, unfortunately).

Notes on off-target cleavage

Off-target effects are a somewhat controversial topic. On the one hand, cleavage at unintended locations in the genome is a major concern, it is important enough to get mentioned in newspaper articles and has led to stockmarket changes and calls for study retractions. On the other hand, if you are working with mice, worms or flies, and you are breeding out for a few generations after your genome edits, then you are probably not worried at all, as unrelated mutations are most likely get removed over time. In this case, you can probably skip this section, just avoid guides with very low specificity scores and possibly use the “same chromosome” filter in column 7 of the guide list to show only off-targets that are on the same chromosome, so less likely to get removed by back-crossing.

For other researchers, e.g. those that use mice without further inbreeding or work on cell cultures or on medical applications, the off-target issue needs to be addressed. Here is a quick summary of what is relatively well-known

about off-targets:

- in very small genomes, like bacteria or viruses, it is unlikely to find similar sequences to a 20mer and off-targets are rarely a problem
- off-target effects primarily depend on the guide sequence. A guide sequence that is entirely different from every 20mer in the genome (e.g. many guides targeting GFP) is extremely unlikely to cut anywhere else than at its target site
- as such, any results you have heard from colleagues about hundreds of off-targets apply to their particular guides and may not be relevant for your guides.
- the more mismatches the off-target has relative to the guide, the less likely it is that the off-target is a problem.
- mismatches that are close to the PAM site make it less likely that the off-target is a problem.
- certain nucleotide changes (e.g. C->T) at certain positions make it less likely that the off-target is a problem.
- off-targets cannot be predicted entirely, these rules are not everything. There still are a handful of known, strong off-targets that are unexplainable. They differ from the guide at many positions, yet are strongly cleaved even in-vitro assays of naked DNA (no chromatin structure)
- a few known off-targets, not more than 3-4 in total, to our knowledge, seem to have one-basepair insertions or deletions relative to the guide
- new high-fidelity versions of SpCas9 can reduce the off-target effects, but not to 0, there still will be some cleavage on some sites (more on this later)

CRISPOR aligns the guides against the genome, reports those locations with a similar sequence in the genome (“predicted off-targets”), ranks these and then based on their scores, ranks the guides.

When predicting off-targets, CRISPOR searches the whole genome (including repeated regions) and allows up to four mismatches. It is important to search repeats, as they are part of the genome, but you may wonder why we chose four mismatches. When we looked at real off-targets, around 90% of them had not more than four mismatches. While you can increase the number of allowed mismatches using the command line version of CRISPOR to five or even six mismatches, the large majority of the predicted locations will be false positives, and most guides would require thousands of PCRs to check all predicted off-targets with more than four mismatches. See Figure 1 of the CRISPOR paper.

CRISPOR then ranks the putative off-target sequences by a score called “CFD”, created by John Doench and colleagues in his 2016 paper. We have evaluated four different scores in (Figure 2 of the CRISPOR paper) and found

the CFD score more accurate than the others.

CRISPOR summarizes all of the putative off-target sites it finds in the genome for one guide into the “guide score”, a formula we copied from the MIT CRISPR website. The “Hsu” score itself was never published, but is based on data from a paper by Patrick Hsu. While we use Patrick Hsu’s formula like the MIT site, CRISPOR searches go up to four mismatches, and the MIT website uses three mismatches and does not search repeats, so the guide scores are not exactly identical. CRISPOR is more sensitive, but in practice most guides with a low score on the MIT website will also have a low score with CRISPOR. The commercial solutions Benchling and DeskGen are also less sensitive, so their scores differ, but ChopChop and Cas-OffFinder should give identical off-target results to CRISPOR.

So how many off-targets should you expect? We studied this question in Figure 3 of the CRISPOR paper. It is based on 30 guides that were tested for off-targets with high-throughput assays by various labs.

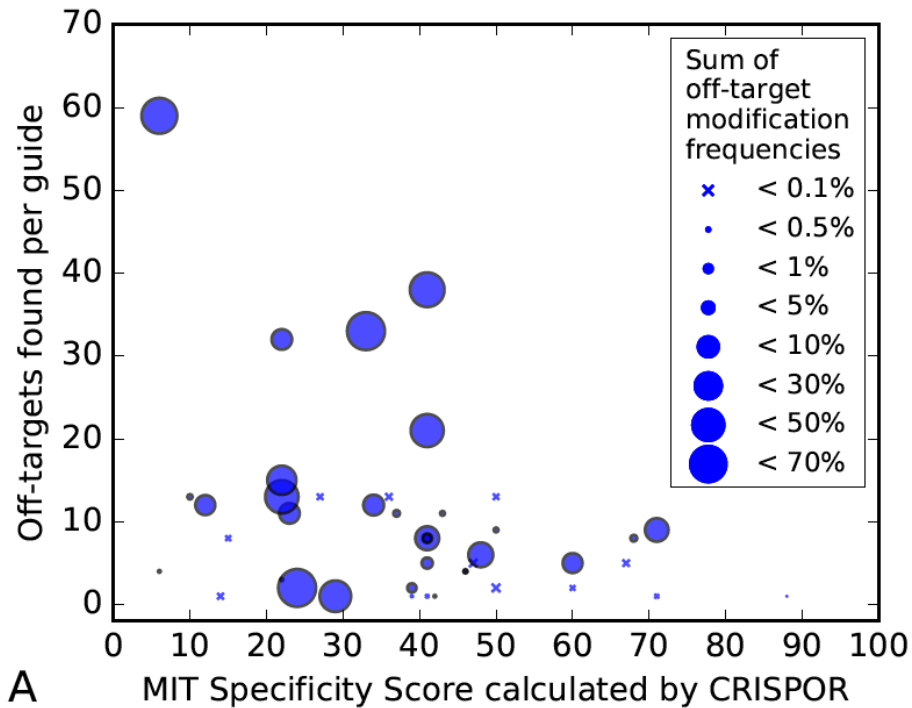


Figure 15: Number and strength of off-targets relative to guide specificity

The figure shows that there are some guides that have virtually no off-targets, but in general, there is a decreasing tendency of strong off-target effects

when the specificity is higher. It also shows that too many tests have been done for low-specificity guides, which may explain why you have heard from colleagues that hundreds of off-target sites are common. It also shows that the accuracy of the prediction is not great, the current predictions have some value, but there is no guarantee that you can avoid off-target effects with a software like CRISPOR alone.

Therefore, if you are worried about off-target effects, you can either test your guide for off-targets with a high-throughput technique in cell cultures (e.g. guideSeq or DiGenomeSeq2), which will require learning a new assay. The alternative is to test predicted off-target sites by a simple PCR and next-generation sequencing (or use TIDE, knowing that it will not find the weaker off-targets). For a highly specific guide, this may involve only a few dozen PCRs and CRISPOR's Off-target primer designer makes their design very easy. This is the other good reason to use only very specific guides: you can screen all predicted off-targets by PCR.

How many off-targets should you test, especially since there are often hundreds of predicted sites and only very few real off-targets among them? There is not enough data to answer this question completely, but here are a few examples (complete histograms):

- the most unspecific guide that we know, VEGFA_site3 from the GuideSeq paper, responsible for the majority of all known off-targets to date, has a MIT specificity score of only 6. GuideSeq found 32 off-targets in Tsai et al, but 4450 off-targets are predicted by CRISPOR. If you sequence only the off-targets with a CFD score over 0.2, then you have to test 274 potential off-targets but still find 24 real off-targets among them.
- the somewhat more specific guide EMX1, MIT score 21, has 11 real off-targets found by GuideSeq and 1026 predicted off-targets. There are only 143 off-targets with a CFD > 0.1, which still catches 10 out of the 11 real off-targets. This is somewhat more than a 96-well plate but still doable.
- the VEGFA guide from the DiGenomeSeq paper has a reasonable specificity of 42. You may not want to test all 719 CRISPOR predictions, but if you limit yourself to the ones with a CFD > 0.1, you have to sequence only 75 and still find 21 out of the 24 known off-targets.
- guides with a higher specificity have very few off-targets and almost all of them will be found by testing predictions with CFD scores > 0.1. This requires testing 30-70 off-targets.

It seems that if your guide has a reasonable specificity score, testing one 96-well plate of primers should find most off-targets.

Should you use a high-fidelity version of SpCas9, eSpCas9 1.1 or HF1? If you

are worried about off-targets, then probably yes. Just do not be surprised that these enzymes are less efficient on many target sites or may not even be able to cleave some target sites at all. If you use them, you may want to try both enzymes or be ready to fall back to the original SpCas9. Kulcsar et al. 2017 have studied this topic in detail and shown that there is a ranking from the best to the worst enzyme in terms of efficiency, but the order was different for every target site they tested.

Notes on on-target cleavage

On-target efficiency, the percentage of cells with a DNA modification of the target sequence, can vary a lot between guides. The general recommendation is to first test several guides, often three, around a target sequence and then use the best one. Quite a few labs have looked at their CRISPR results post-hoc to determine features of their most efficient guides. Initially, these were rules, derived from smaller worm or fly datasets, and lead to guidelines like “prefer guides with Gs at the 3’ end”. Later, pools of thousands of guides were tested in cell cultures or, one by one, expressed in-vitro and injected into zebrafish. The results of these assays were run through machine learning algorithms to build models that predict the efficiency from on the guide sequence alone.

In our CRISPOR paper we conducted the most comprehensive comparison of these prediction models at the time. No matter the model, their value is limited, they can explain only ~40% of the variation between different guides. The unexplained part does not seem to depend on the chromatin structure or the other possible factors that the studies have tested. Therefore, when looking at your own results, if you compare only a few guides, you may not see a significant correlation of your results against the prediction scores. In addition, if you use a very inaccurate method to quantify your efficiency, like the T7 endo assay, many differences will not be visible. Finally, according to a Synthego slide deck, synthetic RNA guides do not seem to have this link between sequence and efficiency at all (watch this space, we will follow up on this soon).

Nevertheless, if you plan to test more than just a few guides, the efficiency prediction scores do enrich for more efficient guides. We found that there is one important parameter when choosing a score, it is the way you express your guides. T7 in-vitro expressed guides - injected into e.g. mouse or Xenopus eggs - behave very differently from U6 guides expressed from plasmids that have been transformed into cell cultures. In our analysis, the score by Moreno-Mateos et al. worked better for T7 datasets, and the one by

Doench et al. 2016 was best for U6 datasets, see Figure 5 of our CRISPOR paper paper.

The influence of the efficiency score on the measured cleavage is shown in the diagram below: the y-axis shows the Doench2016 score range and next to each range, a histogram illustrates how many of the tested guides in this range were among the top 25% efficient guides (green) and how many were among the bottom 25% guides (red) in the whole dataset. Light-red and yellow are the intermediate quartiles. You can see that when guides were tested in the 90-100 score range, very few of them were in the least efficient quartile. In the 10-20 score range, a majority of guides were not very efficient. This is not due to “over-fitting”, because the diagram shows the Hart 2016 experiment, which was not published yet when Doench 2016 was accepted.

The lower diagram shows the same for zebrafish data, but since the dataset from the same authors is plotted here, the results are probably too good and should be somewhat worse in practice, more like the cell culture diagram.

So while you may have heard that efficiency predictions are of little value, we still recommend them. They should save screening time in the long run as they decrease your chances of stumbling over a very inefficient guide. Just make sure that you look at the right score for your expression system.

Frequently Asked Questions

- Why is the MIT score (aka Hsu score) displayed by CRISPOR different from the ones displayed by the MIT website, Desktop Genetics or Benchling?

Because these other tools do not find as many off-targets as CRISPOR. As we have shown in our paper, CRISPOR finds all off-targets with 4 mismatches. The MIT websites does not find all off-targets, because it filters for repeats. Benchling misses some off-targets, not as many as the MIT website but still a high number. We do not know why, but it may also be related to repeats. Desktop Genetics only can find off-targets with up to three mismatches, so their MIT score is also slightly different.

- Why is there a CFD off-target score but not a CFD guide score?

Because the CFD score by Doench et al. is not defined on the level of guides. It is only useful for ranking off-targets. This is why we and others use the MIT score to rank guides.

- Can I score my existing single guides with CRISPOR ?

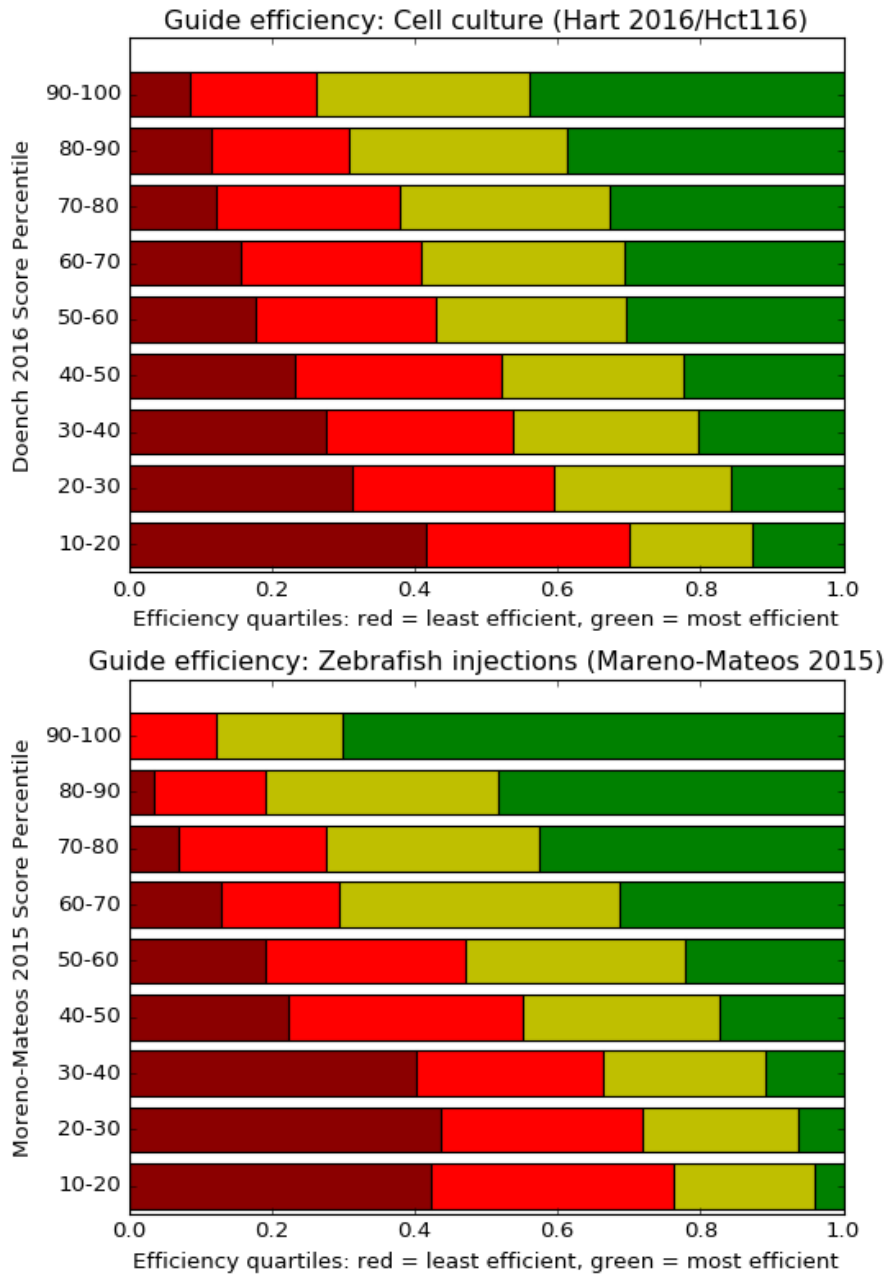


Figure 16: Distribution of efficiency by prediction score bin

Of course, just paste the target sequence, i.e. the guide + the PAM sequence, into the sequence input box.

- Can you add my genome?

Yes, please send us the UCSC genome name or the Ensembl name or the NCBI RefSeq Assembly ID, which starts with GCF_ or GCA_ . We can also add individual FASTA and GFF files, which you can send by URL or as a Dropbox link. We have received more than 150 genomes from users by email until now.

- Can you add my genome and not share it?

Yes, we can add private genomes. People can use them on the website, but they will not be able to download the genome. Just tell us that the genome is private and shall not be downloadable.

- Can you add my genome? By the way, it is 10GB big.

This is tricky. The problem is usually not the size of the genome, but rather the insane number of contigs (often several millions) which makes sequence searches very slow. Searching such a number of sequences takes forever and will make the website unusable for other users. One solution is to reduce the number of contigs, by concatting all smaller ones into a single “chrU” contig, like in the older human genome versions. Another solution is for you to setup your own personal CRISPOR website with just your genome. See the Downloads page at the bottom of the CRISPOR start page to learn more and contact us for help.

- Can I search for 19bp long guides?

Yes, but only with the CRISPOR command line version. It has the `-shortGuides` option. Look for the command line version on our Downloads page. If you have trouble, email us. Because of their lower efficiency, we do not recommend 19bp long guides and do not offer this option on the website. If you think this is a mistake, let us know.

- Does CRISPOR support nickases?

Nickases are a variant of Cas9 that require two guides with targets very close to each other and both have to cut for a successful genome edit. We have received numerous reports that two guides are a lot less efficient than single guides and less specific than the new high-specificity Cas9 enzymes. We have no plans anymore to support explicit two-guide nickase designs. If you think that this is a mistake, let us know, but for now, we do not recommend the double-guide nickase strategy anymore. We rather recommend eSpCas9 1.1 or SpCas9-HF1 and single guides.